

## spirobank™ II Spirometer



### PATIENT MONITORING

#### Unrivalled features in a single device

The **spirobank™ II** is a handheld multifunction spirometer with graphic display. It can be used alone or with the winspiroPRO software connected to a PC. The spirobank™ II does not require calibration and has an internal temperature sensor for automatic BTPS conversion. An internal 6,000 patient memory with Bluetooth® and telemedicine communication features is standard. An optional built-in pulse oximeter makes the spirobank™ II the most versatile and user-friendly handheld spirometer on the market. The oximetry option allows the user to perform oximetry tests such as oxygen qualification, overnight sleep screenings, and real-time oximetry readings. This information can be downloaded to the winspiroPRO software for report generation.



#### Disposable Turbine – Industry First!

The world's first disposable turbine offers outstanding patient sensitivity while providing a reduction in per patient test costs. The disposable turbine reduces the chance of cross-contamination, thus ensuring each patient's safety.



## spirobank™ II Spirometer

- Spirometer with 6,000 spirometry test memory capacity
  - High resolution/high visibility graphic display
  - Patient identification by name or ID code
  - Automatically records the ATS best test with 3 trials
  - Pre-post bronchodilator comparison
  - Spirometry interpretation with test quality assurance
- **Optional Oximeter**
  - Short or long term (overnight) SpO<sub>2</sub> and Pulse Rate measurements
  - Sleep screening with desaturation events recording
  - Walk Test (6MWT) with SpO<sub>2</sub> and Pulse Rate analysis
  - Oximeter option with 1,000 hours (40 days) oximetry recording

## winspiroPRO Software

- On-line PC connection with icon interface
- Real time Flow/Volume and Volume/time curves
- Bronchial challenge with FEV1 dose-response
- Pediatric incentive animations
- Lung age estimation
- Integration with EMR systems
- Data and graphs export via e-mail



FVC Analysis




MVV Analysis



Post Analysis



SpO<sub>2</sub> Analysis

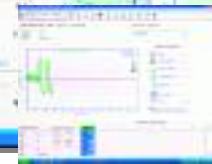

VC Analysis

## Technical Specifications

- **Spirometry Features and Measured Parameters**  
FVC, FEV1, FEV1/FVC%, FEV3, FEV6, FEV1/FEV6%, PEF, FEF25%, FEF50%, FEF75%, FEF25-75%, FET, Vext, FIVC, FIV1, FIV1/FIVC%, PIF, \*FVC, \*FEV1, \*PEF, VC, IVC, IC, ERV, FEV1/VC%, FEV3/FVC%, VT, VE, Rf, ti, te, ti/t-tot, VT/ti, MVV (\* = Best value)
- **Temperature Sensor**  
Semiconductor 32-113°F (0-45°C)
- **Flow Sensor**  
Bi-directional digital
- **Volume Accuracy**  
± 3% or 50 mL/s
- **Flow Accuracy**  
± 5% or 200 mL/s
- **Dynamic Resistance at 12 L/s**  
<0.5 cmH<sub>2</sub>O/L/s
- **Display**  
STN graphic, 128 x 64 pixel
- **Keyboard**  
Membrane, 6 Keys
- **Connectivity and Data Transmissions**  
USB, Bluetooth®, wireless internal modem for acoustic coupling to telephone (telemedicine applications), RS 232 (cable optional)
- **Power Supply**  
4 x 1.5 V, AAA battery
- **Weight**  
6.35 ounces (180 grams) with battery
- **Dimensions**  
Length: 5.8 inches (145 mm)  
Width: 2.4 inches (60 mm)  
Height: 1.2 inches (30 mm)
- **Standards**  
CE-0476, ATS-ERS-Standards, ISO-9001-2000, ISO-13485
- **Oximeter (optional)**
- **SpO<sub>2</sub> Range**  
0-99%
- **SpO<sub>2</sub> Accuracy**  
± 2% between 70-99% SpO<sub>2</sub>
- **Pulse Rate Range**  
30-254 BPM
- **Pulse Rate Accuracy**  
± 2 BPM or 2%

Information subject to change without notice. Please read the instructions for use supplied with the product for detailed instructions, warnings and cautions.

Q-Line to U.S.A. customers call 1-800-348-1842 for emergency and clinical application information 24 hours a day, 365 days a year.

For more information, please call Smiths Medical PM, Inc., at 262-542-3100 or 800-558-2345, or your Smiths Medical distributor.

### Smiths Medical PM, Inc.

Patient Monitoring and Ventilation

N7 W22025 Johnson Drive, Waukesha, Wisconsin 53186 USA

Phone: 262-542-3100 Fax: 262-542-0718 Toll-Free: 800-558-2345

www.smiths-medical.com

Smiths Medical – a part of Smiths Group plc

The Smiths Medical Design Mark and BCI are registered trademarks of the Smiths Medical family of companies.

The symbol © indicates the trademark is registered in the U.S. Patent and Trademark Office and certain other countries.

All other names and marks mentioned are the trade names, trademarks or service marks of their respective owners.

©2007 Smiths Medical family of companies. All rights reserved. spiro II ss rev 00 03/07

smiths medical