

SPIROVIT SP-1

Simply Breathtaking

Spirometry made simple and affordable!

The SPIROVIT SP-1 is a unique solution to simplify complex lung function testing. The high resolution patient incentive screen displays real-time flow/volume loops and detailed measurement matrices. The data can be either printed within seconds or transmitted to PC.

Disposable Sensor SP-150

With SCHILLER's unique SP-150 disposable sensors, the risk of cross-contamination is reduced to its minimum. The small, lightweight sensors are low cost and completely disposable. Simply remove, discard and replace the sensor after each patient performance.

No need for sterilization, easy handling and saving time are maximum benefits for both patient and physician!

Reusable Sensor SP-20

SCHILLER's SP-20 reusable sensors are small, lightweight and very easy to clean. The low cost filters provide the most economical method of performing pulmonary function testing.

Features

- FVC, SVC, MVV, MV and pre/post medication tests with choice of predicted normals
- Simple operation
- Internal memory for up to 100 tests (optional)
- Built-in mains and battery operation
- Ready-to-file reports within seconds on 3.5" (90 mm) paper
- Comparison of predicted and actual values
- Automatic interpretation
- Data transmission via built-in RS-232 interface to PC-based data management program SEMA for longterm storage, retrieval and trending
- Everything in one box ... memory, battery, printer, patient incentive screen
- Affordable, fast, reliable, CE labelled



SCHILLER

www.schiller.ch

The Art of Diagnostics

Answer fax

Yes – I am interested in the innovative products from SCHILLER. Please send me more information about:

- | | |
|---|--|
| ☐ 3/6/12-channel-ECG | ☐ Long-term blood measurement |
| ☐ PC ECG | ☐ Patient monitoring |
| ☐ Ergometry | ☐ Data management |
| ☐ Spirometry | ☐ Telemedicine |
| ☐ Ergospirometry | ☐ |

Name
Clinic, Practice
Function
Address
Telephone
Fax
E-Mail

Asia

SCHILLER Asia-Pacific

10 Jalan SS 3/33
Taman Universiti
47300 Petaling Jaya
Selangor, Malaysia
Phone +603 78 77 5336
Fax +603 78 77 5744
e-mail: schiller@top.net.my

Austria

SCHILLER Handelsgesellschaft mbH

Kampmüllerweg 24
A-4044 Linz
Phone +43 732 70 99-0
Fax +43 732 757 000
e-mail: sales@schiller.at

C. I. S.

SCHILLER Representative Office

Pokrovsky Blvd 3
Office 2202
P.O. Box 27
109028 Moscow, Russia
Phone +70 95 956 29 10
Fax +70 95 956 18 63
e-mail: service@schiller.ru

France

SCHILLER France S.A.

Zac des Luats
58, Route de Champigny
F-94350 Villiers sur Marne
Phone +331 49 41 24 40
Fax +331 49 41 24 49
e-mail: info@schiller.fr

Germany

SCHILLER Medizintechnik GmbH

Rudolf-Diesel-Strasse 14
D-85521 Ottobrunn
Phone +49 89 62 99 81-0
Fax +49 89 609 50 90
e-mail: info@schillermed.de

India

SCHILLER Healthcare India Pvt. Ltd.

D.C.Silk Mills Compound,
'A' Wing, 1st floor
5, Chunawala Estate, Kondivitta Lane
Andheri-Kurla Road, Andheri
Mumbai – 400 059, India
Phone +91 22 826 35 20
Fax +91 22 826 35 25
e-mail: schiller@vsnl.com

Spain

SCHILLER España, S.A.

Avenida Baviera, 12-14 Posterior
E-28028 Madrid
Phone +34 91 713 01 76
Fax +34 91 355 79 33
e-mail: schiller@schiller-es.com

Switzerland

SCHILLER-Reomed AG

Riedstrasse 14
CH-8953 Dietikon
Phone +41 1 744 30 00
Fax +41 1 740 37 10
e-mail: sales@schiller-reomed.ch

Turkey

SCHILLER TÜRKİYE

Bayar Cadd. Paker Sitesi 97B
TR-81651 Kozyatagi-Istanbul, Türkiye
Phone +90 216 372 3730
Fax +90 216 463 7639
e-mail: dtctr@superonline.com


USA

SCHILLER America Inc.

11300 NW 41st Street
Miami, Florida 33178, USA
Phone +1 786 845 06 20
Fax +1 786 845 06 30
e-mail: asjobson@aol.com

Technical Data:

Dimensions: 290x210x69 mm, app. 2.9 kg; 11.4 x 8.3 x 2.7 in, app. 6.3 lbs
Power Supply Requirements: 220-240V, nominal 230V, 50/60 Hz, or 100-115V, nominal 115V, 50/60 Hz mains independent operation with built-in rechargeable battery and charger unit, LED indication for line or battery operation
Battery Charging Time: 15 hours from a completely discharged battery (<3 hours to 60% capacity)
Battery Capacity: 6 hours of normal use
Power Consumption: max. 28 VA
Storage: Memory for up to 100 tests
Liquid Crystal Display: - Liquid crystal display for graphic and alphanumeric representation - Resolution: 192 x 64 dots - Contrast is manually variable via keyboard
Chart Paper: Thermoreactive, Z-folded, width 90 mm/3.6 in, perforation each page
Printing Process: High-resolution thermal printhead, 8 dots/mm
Paper Speed: 25 mm/s
Safety Standards: RL 93/42/EEC, Annex II; EN 60601-1: 1990; pr EN 1441: 1994; IEC 601-1-2
Safety Class: BF according to IEC 601-1
Protection Class: - I according to IEC 601-1, CSA 22.2 and UL 2601 (units with internal battery) - IIa according to RL 93/42/EEC (medical guideline)
Conformity: CE-0124 according 93/42/EEC
Control Panel and Keyboard: User-friendly Rubber Keys, LED indicators
Environmental Conditions: - Temperature, operating: 10° to 40° C/50° to 104° F - Temperature, storage: -10° to 50° C/14° to 122° F - Relative Humidity: 25 to 95 % (not condensed) - Pressure, operating: 700-1060 hPa
Standard Interfaces: 1 x RS-232 serial interface
SPIROVIT SP-1 Basic Unit: Performs four inspiratory and expiratory pulmonary function tests with choice of predicted normals
Accessories: - 1 pneumotach sensor (choice of SP-150 disposable or SP-20 reusable sensor) - 1 noseclip - 1 pack of disposable mouthpieces (for SP-150) or filters (for SP-20) - 1 power cable - 1 pack of paper - 1 operating manual
Hardware Options: Transport case, calibration syringe
Software Options: - SEMA PC-Software for validation and archiving of test data on PC - Memory for up to 100 tests
Measured Values: FVC: FVC, FEV _{0.5} , FEV _{1.0} , FEV _{3.0} , FEV _{0.5} /FVC, FEV _{1.0} /FVC, FEV _{3.0} /FVC, FEF _{0.2-1.2} , FEF ₂₅₋₇₅ , FEF ₇₅₋₈₅ , PEF, MEF ₇₅ , MEF ₅₀ , MEF ₂₅ , FIVC, FIV _{1.0} , FIV _{1.0} /FVC, FIV _{1.0} /FVC, PIF, MIF ₅₀
SVC: SVC, ERV, IRV, TV
MV: MV, RR, TV
MVV: MVV, RR, TV
Prediction Equation: - Adults: ECCS, Crapo, Morris, Knudson, Knudson 76, Austria, Berglund, Finnish, Indian, Composite, Polgar - Children: Quanier & Tammeling, Austria, India, Knudson, Knudson 76, Polgar
- Comparison pre/post medication possible
- Extrapolated predicted values
Presentation: - Flow volume loop - Flow time plot - Measurement values table
Standards Compliance: ATS, ATS 94, OSHA, NIOSH
SPIROVIT SP-150 Pneumotach Flowsensor for Pulmonary Function Testing with Disposable Sensors: Dimensions: 118x36x28 mm, app. 120 g; 4.6x1.4x1.1 in., app. 0.26 lbs
Measuring Method: Pneumotachometer
Measuring Ranges: - Flow: 0 to ±14 l/s - Volume: 0 to ±11 l
Measuring Accuracy: According ATS 1994
Flow Impedance: <0.2 mbar *s/l at 12 l/s
SPIROVIT SP-20 Pneumotach Flowsensor for Pulmonary Function Testing with Reusable Sensor: Dimensions: 125x36x28 mm, app. 160 g; 4.9x1.4x1.1 in., app. 0.34 lbs
Measuring Method, Measuring Ranges, Measuring Accuracy, Flow Impedance same as SP-150
Technical Data are subject to change without notice.

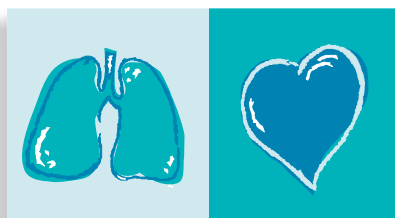

The Art of Diagnostics

SCHILLER

www.schiller.ch